



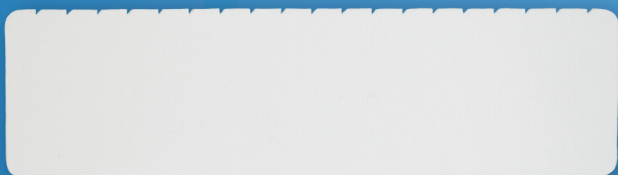


The TZ Medical Adjustable Tensioning System (ATS) is an adjustable aid designed to help facilitate labor induction by creating and maintaining variable tension once the intracervical catheter is placed. The ATS utilizes a soft adjustable tensioning tab that connects directly to the end of the balloon catheter, and a micro hook & loop adhesive pad that is placed along the patient's leg.

Once the balloon is inflated in the cervix, the tensioning tab can be secured to the adhesive anchor pad on the patient, creating lateral tension to assist mechanical dilation of the cervix. Adjusting tension is as simple as pulling directly up on the tensioning tab, and resecuring it to the anchor pad.

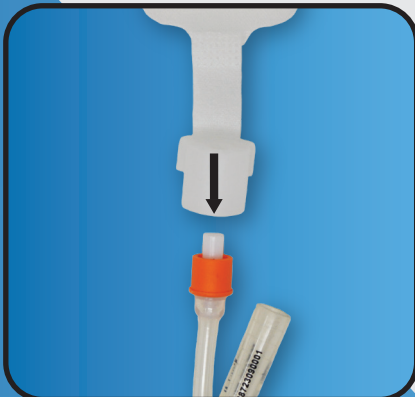

Features

- Provides consistent tension on the balloon catheter to help facilitate cervical dilation
- The adjustable tensioning tab connects directly to the hub of the balloon catheter
- Completely adjustable system to maintain consistent balloon tension through the cervix
- Adhesive pad is skin-friendly and latex-free
- Quick and simple setup

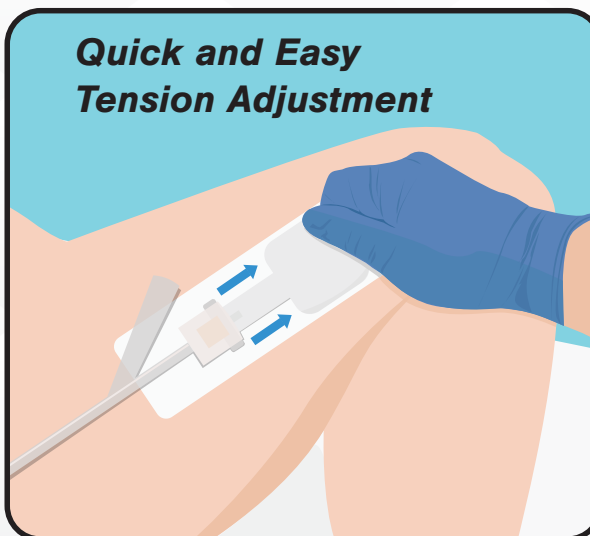

Item Number	Description	UOM
ATS-1005	Adjustable Tensioning System	Box of 10

