



TZ Medical
Sparked by your ideas

PEDIATRIC DEFIB PADS




TZ Medical's Pediatric Defibrillation Electrodes boast the same clinical benefits as our adult pads, however, are designed for your adolescent patients. Our Pediatric Defibrillation Electrodes are constructed to adhere and conform to the patient's body.

Our multi-function electrodes are designed for EP and Cath Labs, Cardiac Surgery, Critical-Care Areas and other related departments. They are designed for defibrillation, synchronized cardioversion, external pacing and ECG monitoring. Compatible with biphasic defibrillators, our electrodes offer direct connections to a variety of defibrillators, such as Zoll® and Medtronic® QUIK-COMBO™. We also offer adapters, allowing multiple brands of defibrillators to use the same electrodes.

Features

- 60" lead wires for maximum convenience and utility
- Hydrogel maintains excellent adhesion to the patient's skin
- Radiotransparent and radiotranslucent options available

Radiotransparency

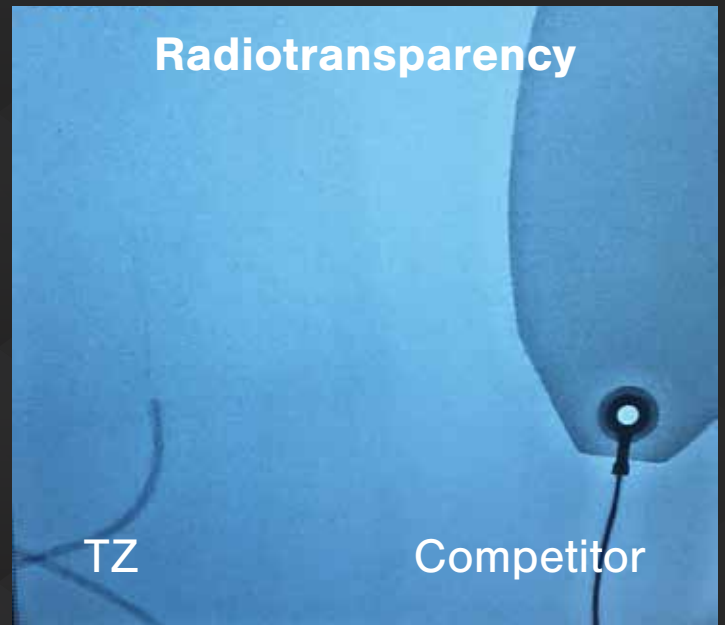

TZ

Competitor

FIND OUT MORE



tzmedical.com



800.944.0187




info@tzmedical.com

19808 R0.10

Pad Placement Guide



TZ Medical
Sparked by your ideas

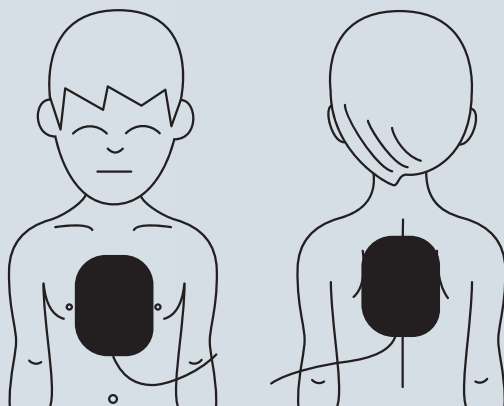

Lead Wire Length

TZ Sterile P-122-T1: 120"

TZ Electrodes: 60"

Competitor Electrodes: 42"

Competitor Electrodes: 24"

Medtronic® M1 Connector



Universal T1 Connector



Zoll® Z1 Connector



Product Description	Size	Medtronic® ¹	Universal ²	Zoll® ³
Pediatric Radiotranslucent pads	C	P-112-T1 PAD 1111 Adapter*	P-112-T1	P-112-T1 PAD 101PD Adapter*
Pediatric Radiotranslucent pads, sterile	C	P-122-T1 PAD 1111 Adapter*	P-122-T1	P-122-T1 PAD 101PD Adapter*
Pediatric Radiotransparent pads and leadwires	C	P-212-M1		P-212-Z1

Pad Dimensions:
C=4.25"x2.876"

***Adapters are reusable**

1 – M1 Medtronic models are comparable to Medtronic® Physio-Control QUIK-COMBO™ electrodes.
 2 – T1 Anderson models allow for multiple brand of defibrillators to use the same electrodes. (adapter needed)
 3 – Z1 Zoll Models are comparable to Zoll® multifunction electrodes.

For Pediatric Use: ≥25kg (55lbs)

FIND OUT MORE



tzmedical.com



800.944.0187



info@tzmedical.com

19808 R0.10