




TZ Medical
Sparked by your ideas

PATIENT QUICK START GUIDE

TRIDENT™ PRO MONITOR


Step 1: Electrodes

Connect color coded snap connectors to the adhesive electrodes.

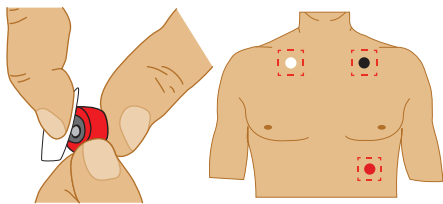


Step 2: Skin Prep

Remove all hair from electrode sites with a razor and scrub the site with warm soapy water. Rinse and dry thoroughly.

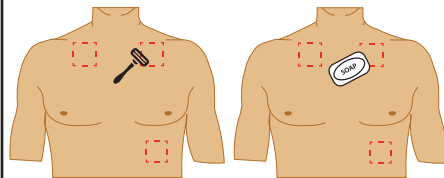


Step 3: Application

Remove backing from the adhesive electrodes and apply to the prepped and dried sites. Make sure to press the border firmly and match the locations shown below.

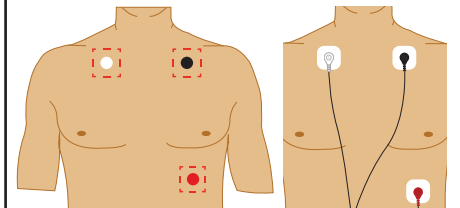

Step 4: Connect Cable

Firmly connect the cable to the device in the orientation shown. Ensure that ECG signals appear briefly on the screen.

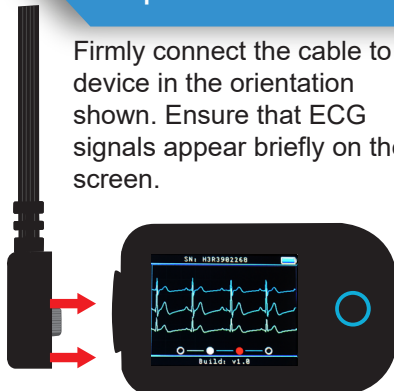


Step 5: Marking Event

To record a symptom, press the Plus button to display the New Event screen. Use the arrows on the bottom left and right of the display to follow the instructions.



Step 6: Record Event

Use the up and down arrows to highlight the symptom and activity you would like to record.

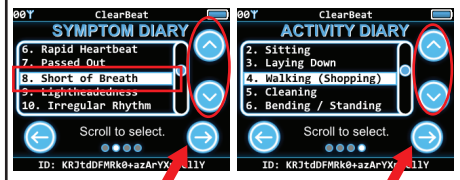


Step 7: Confirm Event

Once the symptom and activity have been selected, press the plus button to confirm the event and continue monitoring.




Step 8: Charging

Charge your device by unplugging the lead cable and plug in the regular charging cable to the device and nearest power outlet until the device is fully charged.

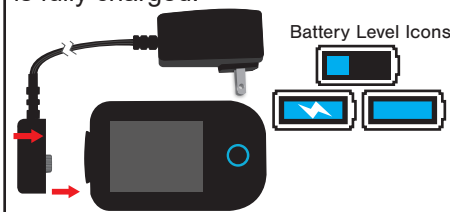

Step 9: End of Study

Power off the device by pressing the power icon in the lower-left of the Done Monitoring screen. Press the power icon to confirm shutdown, and return your monitoring device as soon as possible.

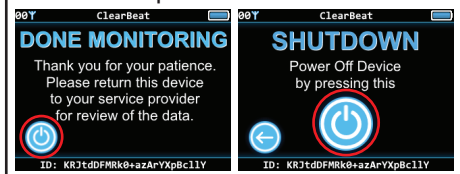

For clinical support, please contact your medical provider.

For technical support related to your cardiac monitoring device, please call **1-888-999-0233** or email **clearbeatsupport@tzmedical.com**. Include your name and preferred contact method (phone, email, other), and a support representative will respond as soon as possible.

TRIDENT™ PRO

What you need to know while monitoring...



Avoid fully submerging or exposing device to excessive amounts of water or moisture. The device is water resistant, not waterproof.



Allow sites to dry completely and remove all residue before applying adhesive electrodes.



During prolonged inactivity, your monitor's display screen will hibernate to conserve battery life. Press the circular, tactile power button to wake the display and access the monitoring screen.



If located in an area with poor cellular service, more frequent charging may be necessary to avoid complete battery depletion.



Make sure to connect all electrodes prior to connecting the cable to the device. If tracings do not display correctly, ensure that the colored snap connectors are securely attached to the electrodes and at the proper locations.



Charging the device at a minimum of every 24 hours is recommended to avoid complete depletion of battery life.

