



**TZ Medical**  
Sparked by your ideas

# STERILE DEFIB PADS




TZ Medical Sterile Defibrillation Electrodes are constructed to adhere and conform to the patient's body, providing superior current distribution, which results in better clinical outcomes.

Our multi-function sterile electrodes are ideal for a variety of procedures, including TAVR, AVR, Mini-Maze, ICD Implants, Laser Lead Extraction and much more. They are designed for defibrillation, synchronized cardioversion, external pacing and ECG monitoring. Compatible with biphasic and monophasic defibrillators, our sterile electrodes offer direct connection to Zoll® defibrillators. We also offer adapters, allowing multiple brands of defibrillators to use the same electrodes.



## Radiotransparency

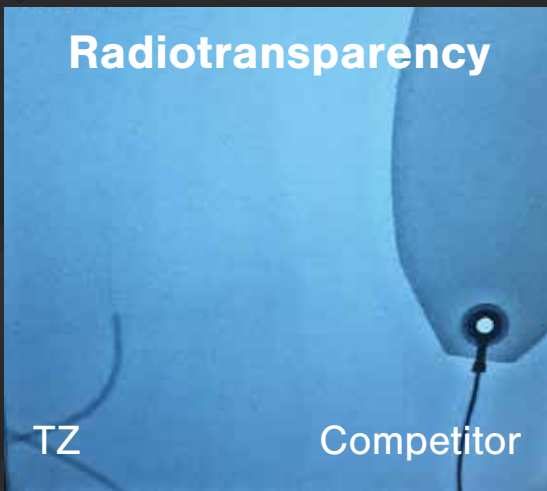

TZ

Competitor

## Features

- Lead wire lengths available in 60" and 120" for maximum convenience and utility
- Hydrogel maintains excellent adhesion to the patient's skin
- Unique edge effect technology disperses energy uniformly, reducing burns
- Radiotransparent and radiotranslucent options available

**FIND OUT MORE**



[tzmedical.com](http://tzmedical.com)



800.944.0187




[info@tzmedical.com](mailto:info@tzmedical.com)

13622 R2.5

# Pad Placement Guide



**TZ Medical**  
Sparked by your ideas

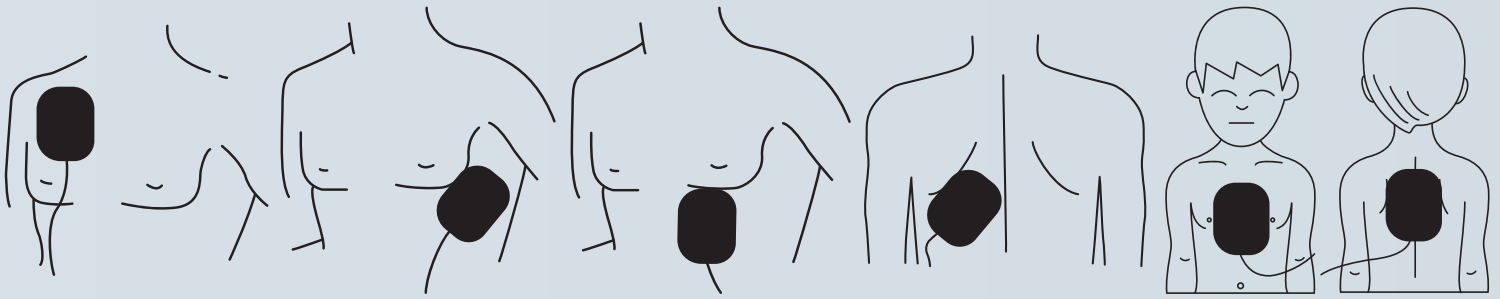

## Lead Wire Length


**TZ Sterile P-121-T1, P-122-T1: 120"**

**TZ Sterile P-224-Z1, P-224-T1: 60"**

**Competitor Electrodes: 24"**

Universal T1 Connector

Zoll® Z1 Connector

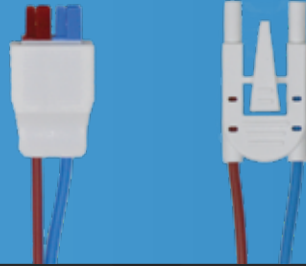

Product Description	Size	Medtronic®	Universal <sup>1</sup>	Zoll® <sup>2</sup>
Adult Radiotranslucent, sterile	A	P-121-T1 PAD 1111 Adapter*	P-121-T1	P-121-T1 PAD 101PD Adapter*
Adult Radiotransparent pads and leadwires, large posterior pad, sterile	A/B	P-224-T1 PAD 1111 Adapter*	P-224-T1	P-224-Z1
Pediatric Radiotranslucent pads, sterile	C	P-122-T1 PAD 1111 Adapter*	P-122-T1	P-122-T1 PAD 101PD Adapter*

**Pad Dimensions:**

**A=6"x 4.25"    B=6"x9"    C=4.25"x2.876"**

**For Pediatric Use:** ≥10kg (22lbs)  
AED or AED Mode: ≥25kg (55lbs) ≥8 years old

**\*Adapters are reusable**

1 – T1 Anderson models allow for multiple brand of defibrillators to use the same electrodes. (adapter needed)

2 – Z1 Zoll Models are comparable to Zoll® multifunction electrodes.

**FIND OUT MORE**



[tzmedical.com](http://tzmedical.com)



800.944.0187



[info@tzmedical.com](mailto:info@tzmedical.com)

13622 R2.5