


# ACCUFLATE INFLATION DEVICE

The TZ Medical AccuFlate is a 20ml disposable device designed to accurately control and monitor pressure during inflation and deflation of dilation balloon catheters. Our patented technology allows the Inflation Device to generate pressure up to 30 atm with accuracy of  $\pm 1$  atm.



**Patented self-locking mechanism for easy one-handed locking control, as well as gradual and rapid pressure release ergonomic handle.**

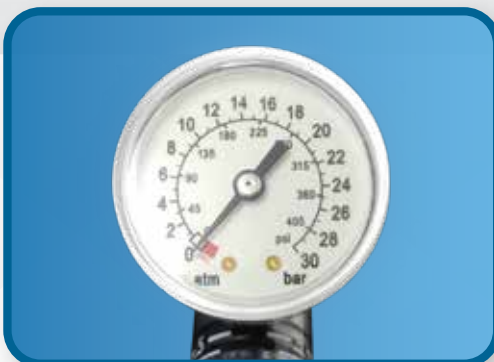

**Precision manometer with a maximum pressure of 30 atm.**



**Optional, male rotating luerlock, 3-way stopcock with "ON" or "OFF" tap handle.**

Item Number	Description	UOM
IDD-0001	Balloon Inflation Device	Box of 5
IDD-0010	AccuFlate with PTCA Kit	Box of 5 sets
IDD-0020	PTCA Accessory Kit with Tuohy-Borst Valve, Guidewire Insertion Tool, & Torque Device	1 each



**PTCA Accessory Kit also available**

**FIND OUT MORE**

Patent: CN103028153A



[tzmedical.com](http://tzmedical.com)



800.944.0187



[info@tzmedical.com](mailto:info@tzmedical.com)

19337 R1.2