

# StandTall™

10 cm, 15 cm, 25 cm



Take advantage of the benefits of Left Radial Artery Access from a comfortable working position. The StandTall™ is an intuitively designed external vascular sheath extender that conveniently redirects workflow to the right side of the procedure table. As a result, you can operate with your preferred dominant hand, maintain optimal posture, and reduce musculoskeletal stress and fatigue.

## Making Left – Right for You

StandTall™ drastically improves both patient and physician positioning for left radial artery procedures performed from the right side of the procedure table. StandTall™ enables the physician to work left radial access cases comfortably, without having to hunch or stoop over, and avoids the hassle of having to revert to femoral access and the various increased risks that going femoral entails. With StandTall™ the physician realizes the benefits immediately through:

### Physician benefits

- Decreased fatigue and musculoskeletal stress
- Improved site access and stability
- Streamlined workflow
- Allows for a more ergonomic posture for the physician from the right side of the table

### Patient benefits

- Greater comfort with left arm positioning
- Quicker recovery and reperfusion times
- Enables convenient access to non-dominant hand (90% of patients are right handed)
- Improved overall procedure satisfaction

### Hospital benefits

- Reduced cost vs. femoral approach
- Reduced post-op patient management
- Fewer physician work-related absences due to MSK stress or injury

DESIGNED TO  
**OPTIMIZE LEFT  
RADIAL ACCESS**



**TZ Medical**  
Sparked by your ideas

## Simple set up, easy to use


StandTall™ can easily be set up in less than 2 minutes.

1. Review patient anatomy and conditions to establish StandTall™ length.
2. Flush StandTall™ and then attach to the vascular access sheath via the universal sheath adapter.
3. Position the adhesive clasp to secure StandTall™ in the optimal position and adhere it to the patient or sterile drape.
4. Shape StandTall™ into the desired position and secure it onto the adhesive clasp of the distal palm.
5. Pronate and position the left arm across the patient to the right side of the table and secure arm to prevent arm drift.

 Without StandTall™



 With StandTall™

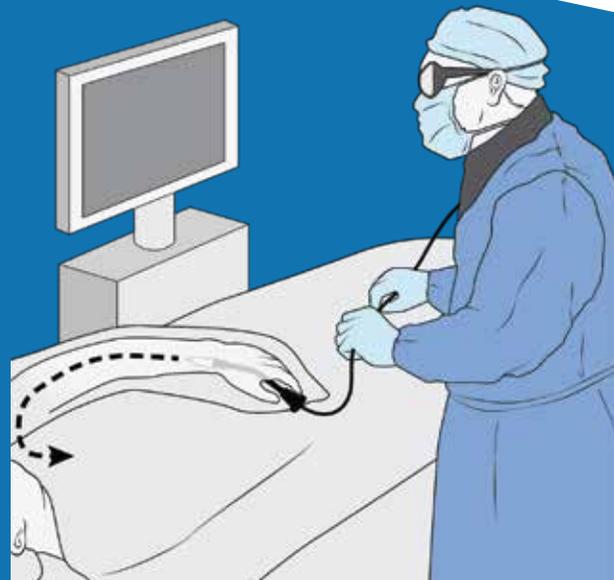

Item #	Description	UOM
RD-ST0010L	External vascular sheath extender 10 cm length	5/box
RD-ST0015L	External vascular sheath extender 15 cm length	5/box
RD-ST0025L	External vascular sheath extender 25 cm length	5/box

