




SUCTION TIPS

TZ Medical's single-use Suction Tips provide optimal control when removing fluid and debris from confined surgical sites. With a wide variety of sizes and lengths, each suction tip is ergonomically designed with a control valve on the handle allowing the operator to regulate the intensity of the suction from a convenient location. The malleability of the aluminum shafts provides better positioning and can be manipulated into multiple different angles. The suction tubes also feature a thin wire or stylette that can be used to remove trapped tissue or debris and provide support for the flexible shaft.

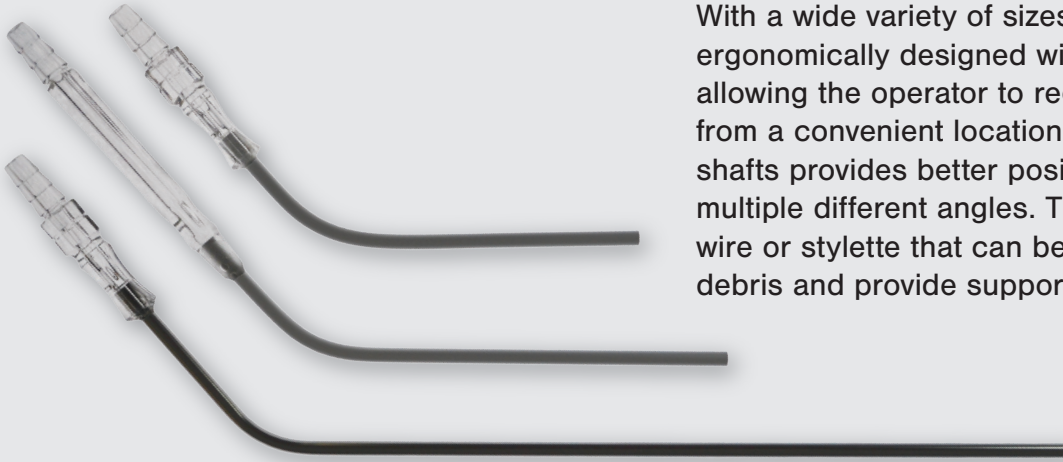

Features

- Single-use and sterile to prevent cross-contamination
- Malleable aluminum shaft provides optimal flexibility and can be molded into any desired position
- All suction tips come with a stylette or thin rod to clear any trapped tissue or debris and prevent obstruction
- Built-in control valve on the handle allows for continuous or intermittent control of suction
- Multiple French sizes and lengths available

Item Number	Description	UOM
STIP-1011-08	Frazier Suction Tip - 6F, 8 cm	Box of 50
STIP-1011-11	Frazier Suction Tip - 6F, 11 cm	Box of 50
STIP-1013-08	Frazier Suction Tip - 8F, 8 cm	Box of 50
STIP-1013-11	Frazier Suction Tip - 8F, 11 cm	Box of 50
STIP-1015-08	Frazier Suction Tip - 10F, 8 cm	Box of 50
STIP-1015-11	Frazier Suction Tip - 10F, 11 cm	Box of 50
STIP-1015-23	Frazier Suction Tip - 10F, 23 cm	Box of 50
STIP-1017-08	Frazier Suction Tip - 12F, 8 cm	Box of 50
STIP-1017-11	Frazier Suction Tip - 12F, 11 cm	Box of 50
STIP-1017-23	Frazier Suction Tip - 12F, 23 cm	Box of 50
STIP-1019-08	Frazier Suction Tip - 14F, 8 cm	Box of 50
STIP-1019-11	Frazier Suction Tip - 14F, 11 cm	Box of 50
STIP-1019-23	Frazier Suction Tip - 14F, 23 cm	Box of 50

